

Lateral flow chromatographic immunoassay | Point-of-care | Easiest handling

Cellex qSARS-COV-2 IgG/IgM Rapid Test

The 15-minute ideal diagnostic test
for SARS-COV-2 Viral infections

- Sample: 10µL Serum/Plasma/Whole blood Specimens.
- Fast and Easy to use: 15 minutes per test, easy to recognize results and no equipment needed
- Affordable for different health systems: Test can be used anywhere, no specialized training is required.
- Accurate and reliable results, robust performance.
- Fast results allow physicians to make immediate evaluation in conjunction with CT scans and other tests.




We Are In The Business
Of Saving Lives



Cellex qSARS-COV-2 IgG/IgM Cassette Rapid Test

3 Simple Steps

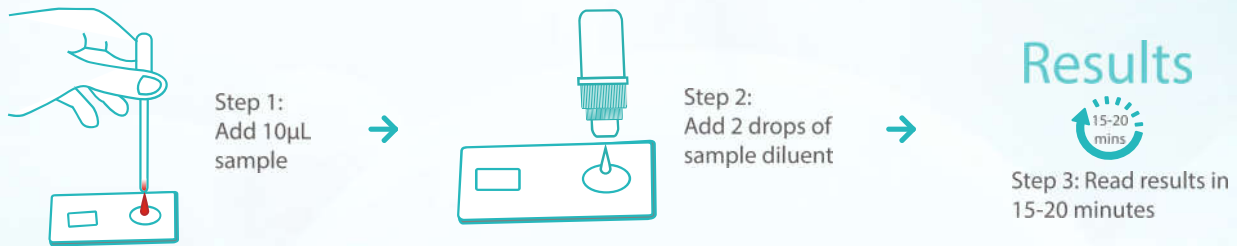

Clinical Performance

A total of 378 samples were collected for this clinical evaluation test, of which 128 confirmed positive samples from COVID-19 makeshift hospital and 250 negative samples collected prior to September 2019

Clinical Test		qSARS CoV-2 IgG/IgM Rapid Test				Sub
		IgG+	IgG-	IgG+	IgG-	
Clinical Status	Pos.	65	46	9	8	128
	Neg.	0	5	4	241	250
Subtotal		65	51	13	249	378

Sensitivity 93.75% (95% CI: 88.06%-97.26%) Specificity 96.40% (95% CI: 92.26%-97.78%)

Interpretation of Result

